


Switch Adapted Water Gun

USER QUICK GUIDE

Introduction

This is a switch adapted water gun that is compatible with 3.5mm switches. This device may be beneficial to someone with limited hand function that is interested in using a water gun.

Features

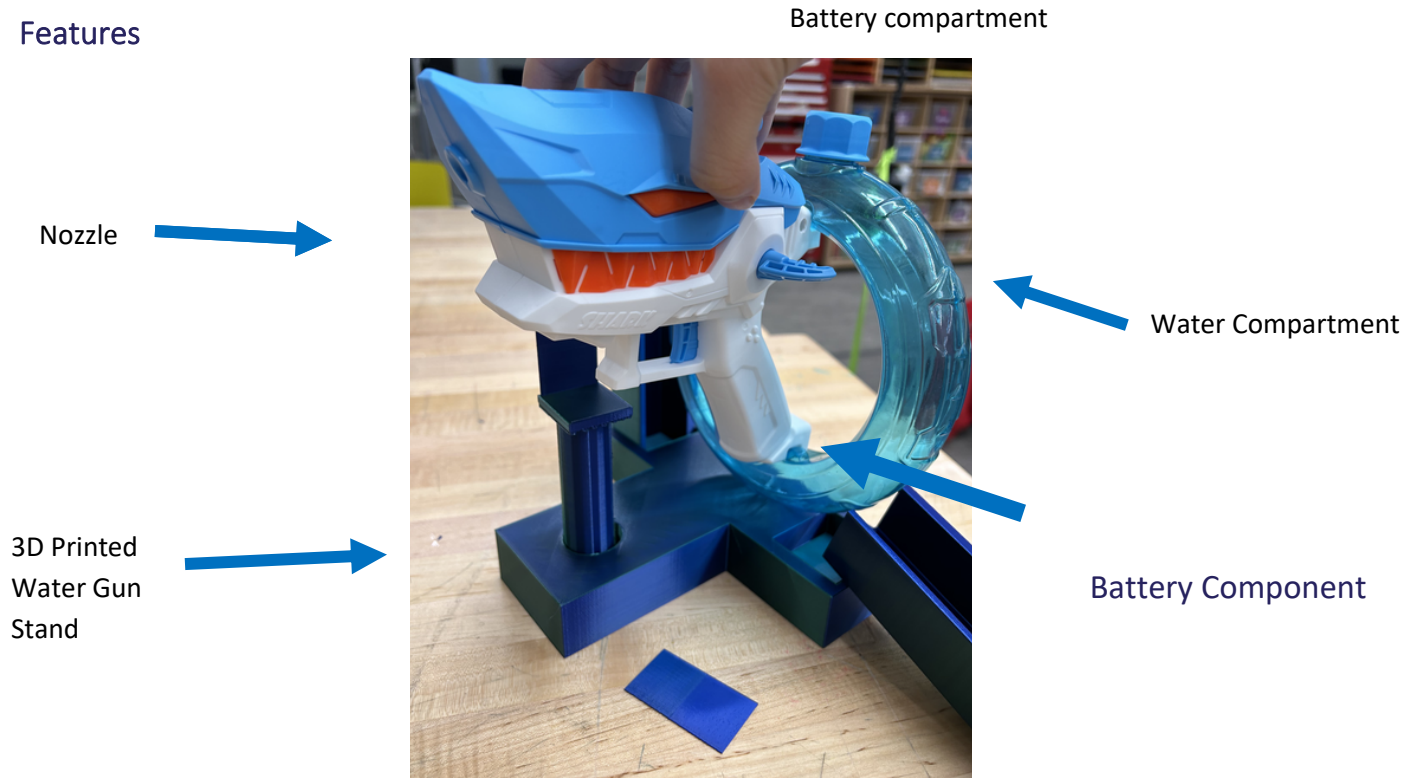

Usage

The Switch Adapted Water Gun can be used by first filling it up with water. Next, plug in a 3.5mm switch into the back of the water gun. Once the switch is pressed, the water gun is activated.

Compatibility

The Switch Adapted Water Gun should be compatible with any 3.5mm switch.

Specifications

Length (mm)	Width (mm)	Height (mm)	Mass (g)
193	87	130	480

*The mass given is the mass of the empty water gun.

Cleaning

The device can be wiped down with a damp cloth. When not in use, empty out the water compartment in the gun.

Switch Adapted Water Gun

USER QUICK GUIDE

Changing the Batteries

The batteries can be accessed by unscrewing all the screws located on the gun's surface. Once they have been removed, open up the gun and access the white cylinder which contains the rechargeable batteries and the charging cable. The batteries can be charged with a usb. The battery shouldn't be charging for a prolonged time.

